

ADDITIONAL INFORMATION

Professional Hydrosana Ionic Detox Foot SPA Machine

A thousand miles begins with a single step



Reflexology can detox, beauty. It is a necessity for everyone!

Product detail:

Name: LCD double belt Detox foot spa machine

Voltage: 110V/220V

Rated power: 60W

Current: 0.5—2.5A

Packaging: 28CM * 42CM * 16CM (L*W*H)

Weight: 6kg

Warranty: One year



Ion detox machine function:

1. Improve the immune system;
2. Improve memory and improve sleep;
3. To improve the reproductive health status;
4. Clean the inside of the body organs;
5. promote liver detoxification;
6. To promote liver, kidney detoxification and discharge parasites;
7. The discharge of heavy metals;
8. relieve pain and stress, stay energized;
9. Improving wrinkles, acne and other skin conditions;
10. reduce body water retention.

The below can be for your reference after many tests:

Yellow-green= detoxifying from the kidney, bladder, urinary tract, female/prostate

Orange= from joints

Brown= liver, tobacco, cellular debris

Black= liver, gallbladder

Dark Green= gallbladder

White foam= mucous from lymph

White speckles= most likely yeast

Black speckles= heavy metals

Red speckles= blood clot material



Product image 1



Product image 2



Product image 3



Product image 4



Regain Health
& Vitality



Heavy Metal Detox Program

Reasons for doing the Heavy Metal Detox Program for YOUR Body!

Hheavy Metals are so pervasive in our society today. There is no longer a question of whether or not we are toxic, but rather what our level of toxicity is.



With everyday living in today's polluted and stressful environment, Heavy Metals can be absorbed from: Industrial Pollution of Air and Water, Lead Pipes, Fumes from Lead-Based Paints, Automobile Exhaust, Cosmetics, Latex Wall Paint, Some brands of Tea, Cheese, Pharmaceutical Drugs, Aspirin, and Table Salt, etc.

Mercury has been found in some vaccines, dental fillings, and fish. Cadmium is found in cigarette smoke, organ meats, and fungicides sprayed on Fruits & Vegies. Aluminum is found in Antiperspirants, Toothpaste, Baby Powder, Cosmetics, Antacids, and Infant Formulas.



The body's excretion pathways may be able to counteract these negative influences for a time, yet become overwhelmed after long term exposure resulting in the symptoms below. This build up of Heavy Metals throughout the body compromises our health & appearance.



Symptoms linked to Heavy Metal toxicity include:



Other symptoms include:

Allergies, Diarrhea, Gastrointestinal Complaints, Parkinson's Disease, Frequent Colds, Numbness, Pain, Skin Problems, Forgetfulness, Tingling, Weight gain or loss, Bloating, Dizziness, Mood Swings, Anxiety, Difficulty waking up, Stress, Bags under the eyes, and Hormone Imbalance.

**Don't let Heavy Metals threaten the Quality of Your Life
Look & Feel Great Inside & Out!
Heavy Metal Detox Successes!**

I have had Chronic Fatigue for 15 years and it was getting worse. After doing a Heavy Metal Detox Program for about 6 months now, I am feeling so much better and more active than I have been in years. I never realized how much effect Heavy Metals could have on my health. Everyone should do this!

60 year old male

I felt miserable, gained weight very easily and had thyroid problems for a long time. It was very hard to lose weight and when I did, it kept coming back. I felt like I was suffering from Depression and was afraid I would have to start taking Anti-depressants. I am so happy now that I did a Heavy Metal Detox program. I am so much happier! I wasn't really depressed after all. I am on a better diet now and not concerned about my weight anymore. Enjoying my life a lot more!

35 year old female

I was always very active, played basketball and ran track. Something changed and I didn't have the stamina I had before. I was getting tired more easily and was afraid I would not be able to participate in sports activities anymore. I had the mercury fillings removed from my mouth a felt a little better, but when I did a Heavy Metal Detox program, that was the icing on the cake. Have much more energy now. So glad that I did this!

20 years old male

Far infrared waist belt use



01. connect the belt with the machine

**02. Press the RED color ON/OFF
switch on the back of the machine
to activate the belt.**



The operating skills

In spite of any countries, there exists the difference of water quality. Only add a few salt, it can promote the ion-cleanse to conduct electricity. The capability of conducting mainly lies on the water electric ability. There exists the capability of conducting more or less in any kinds of waters except distilled water. The capability of conducting of the water determines the ion effects and impacts the function of detoxification. But for maintaining the machine life, please avoid adding salt.

The most normal volume display is 15 V, 0.5 to 2.5 ampere. Don't add too much salt, be noted that if the number of "A" display shows more than 2.8, then the machine will shut automatically.



Five minutes later

Ten minutes later



Fifteen minutes later

Twenty minutes later

Thirty minutes later

Why do we need to preserve our health?




Note



Dangerous

Why do we need to preserve one's health

A global survey from World Health Organization show us that there are only 5% in the world who truly healthy, patients account for about 20%, the remaining 75% belongs to the "sub-health" crowd.

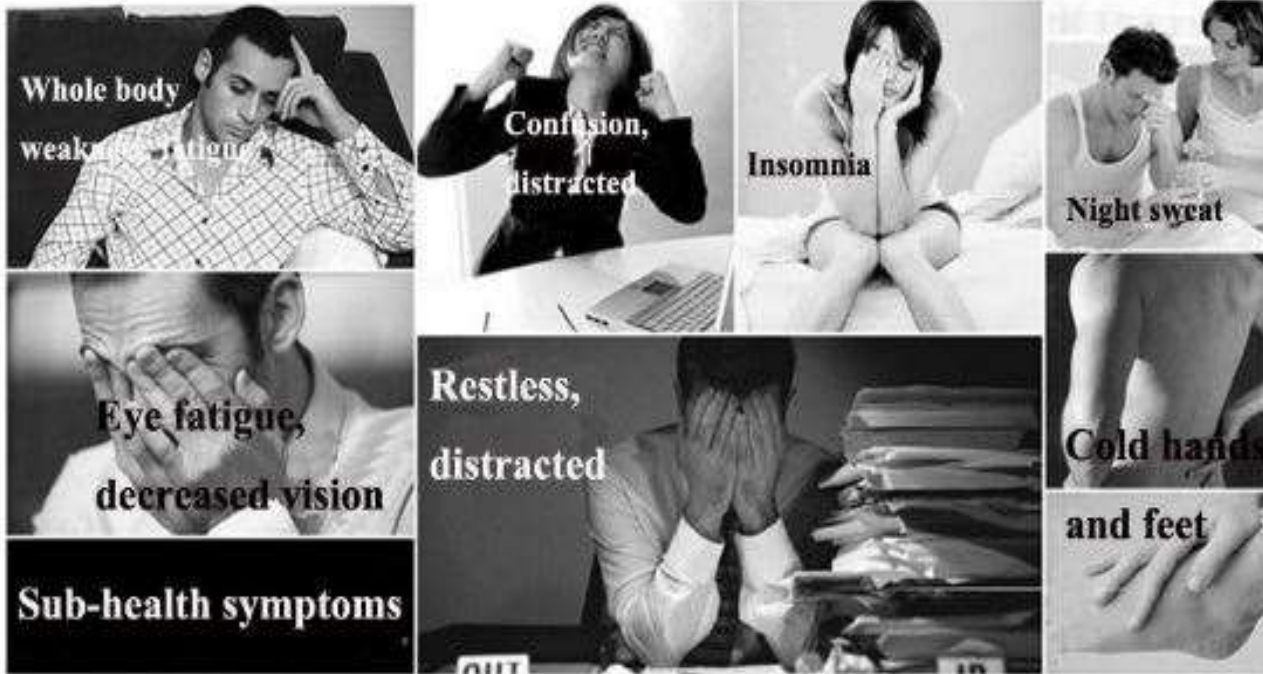

Blood viscosity detection (red blood cell aggregation and dynamic)





Product Accessories:

Array:1

Power source:1

Wrist Belt:1

Aluminum Box:1

Far Infrared Belt:1

Instruction:1

Packaging:

Main machine	Aluminum case	Colorful carton	Array	Wristband	Far infrared waist belt:
					
1pc	1pc	1pc	1pc	1pc	1pc
TENS	Line	Power cable	Certificate of conformity	Warranty card	Manual
					
1pc	1pc	1pc	1pc	1pc	1pc



FOR MORE INFORMATION TO INCLUDE:

1. OTHER DOCUMENTATION

Please email us at: bizdev@freeworldimports.com